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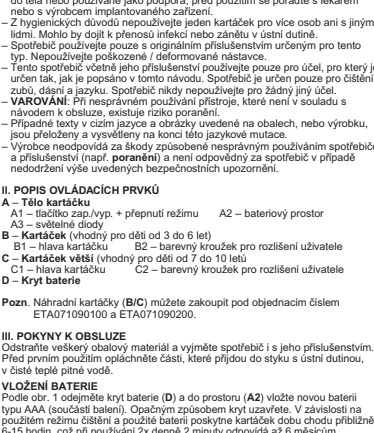
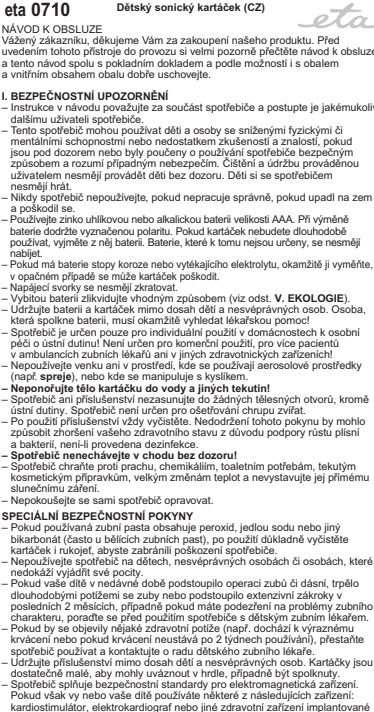
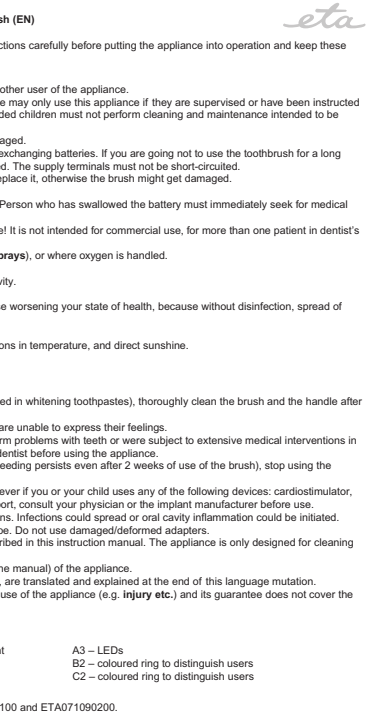
[illegible]

which come in contact with the oral cavity must be rinsed in warm drinking water.

the compartment **(A2)**. Close the cover in the opposite way. Depending on the number of hours of operation which corresponds to about 6 months when used twice a day.

the brush, proceed the other way round. Do not switch on the appliance immediately. Put the toothbrush to the child's teeth at a slight angle towards the front. When you press the button once, you will start the white light mode, catching the tooth to the regular cleaning mode and it will activate colorful lights for more than 2 minutes. When the appliance is switched on, move the brush slowly along the teeth in the mouth. Clean each tooth for about 2-3 seconds. The brush operating time is 2 minutes. If you wish to continue cleaning, simply switch on the brush again (Fig. 2B); the flashing indicator light on the brush and a short interruption of the sound will be heard.

(B1/C1) and the brush body **(A)** under running water and wipe the brush body with a wet cloth for cleaning! Allow the brush body **(A)** to dry properly before its use out of reach of children and incapacitated persons.

the packaging, components and accessories as well as recycling symbols. The appliance may not be disposed off together with communal waste. Hand disposal of the product will help preserving valuable natural resources and avoid waste from improper waste disposal. For more details, please contact the manufacturer or your local waste management authority. Remove battery from the device to avoid the danger of fire!

Accessories for the appliance must be used according to the instructions for use. For repair or replacement of the appliance, please contact the manufacturer or your local service center.

www.ata.cz/servis.eu.