

Ústní sprcha • NÁVOD K OBSLUZE	(CZ)	3-9
Ústna sprcha • NÁVOD NA OBSLUHU	(SK)	10-16
Oral Irrigator • INSTRUCTIONS FOR USE	(EN)	17-23
Szájzuhany • HASZNÁLATI UTASÍTÁS	(HU)	24-30
Irygator dentystyczny • INSTRUKCJA OBSŁUGI	(PL)	31-37
Munddusche • BEDIENUNGSANLEITUNG	(DE)	38-43

AquaCare II

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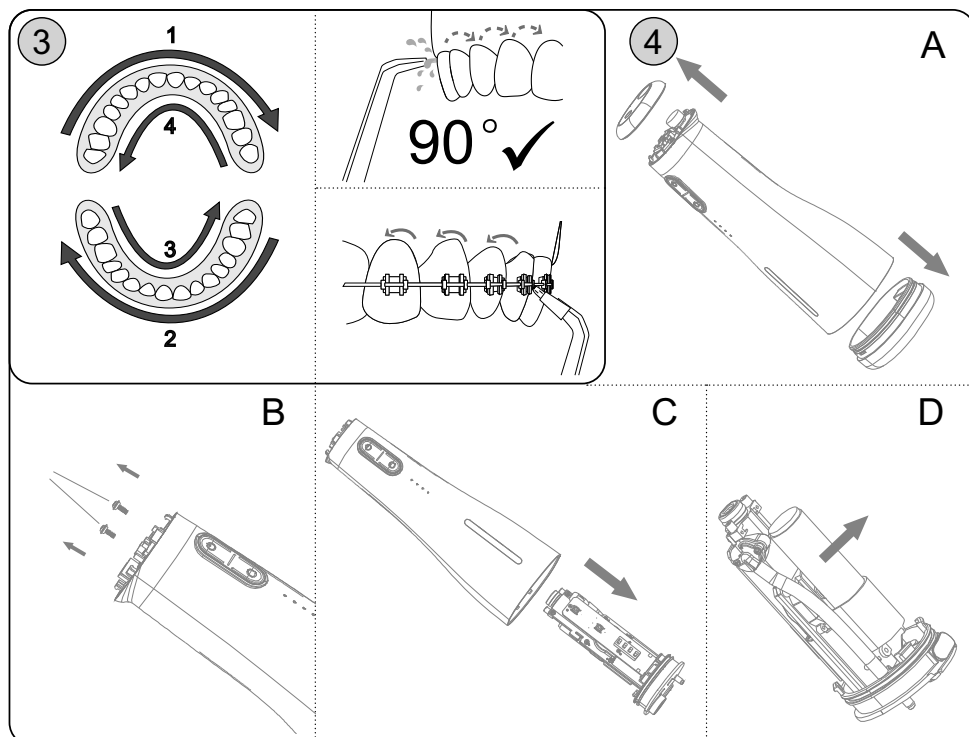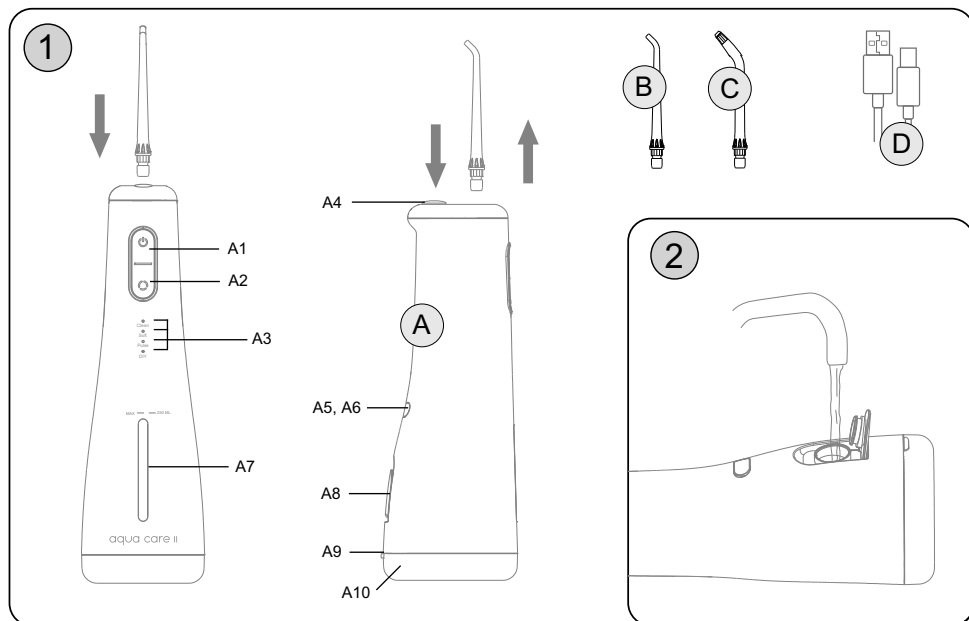Dear customer, thank you for purchasing our product. Please read the operating instructions carefully before putting the appliance into operation and keep these instructions including the receipt and, if possible, the box with the internal packing. This appliance is designed to mechanically remove dental plaque. It is easy to use and enhances your dental hygiene. The oral irrigator features cutting-edge technology in an efficient and ergonomic appliance which is an ideal tool for your dental care.

I. SAFETY PRECAUTIONS

- Instructions in the Instructions Manual are part of the appliance and shall be provided to all users of the appliance.
- Please check if the voltage specified on the rating plate complies with the voltage in your electrical socket. The charging base must be connected to electric socket in the electrical mains complying with the applicable standards.
- This appliance may be used by children aged 8 years and older and by persons with reduced physical, sensory or mental abilities or lack of experience and knowledge only if they are supervised or have been instructed about safe use of the appliance and understand potential risks involved. Children must not play with the appliance. Cleaning and maintenance carried out by the user must not be carried out by children. Children under 8 years of age must keep away from the appliance and its power cord.
- Before replacing accessories or accessible parts that are moving during operation, before assembly and disassembly, before cleaning or maintenance work, make sure to switch OFF the appliance and unplug it from the mains!
- Never use the appliance if its power cord or adapter is damaged, if it does not work properly, if it fell down and was damaged, or if it fell into water. In this case take the appliance to a specialised service centre to check its safety and proper functionality.
- The appliance and the connected adapter must not be used at places where it could fall into bath tub, wash basin or swimming pool. If the appliance fell into water anyway, do not take it out! First of all, unplug the adapter from the socket and only then take the appliance out. Then bring the appliance to a specialised service centre to check its safety and proper functionality.
- If the adapter's power cord is corrupted, it has to be replaced by the manufacturer, manufacturer's service technician or similarly qualified person to prevent hazardous situations.
- When leaving the appliance unattended, before filling the reservoir with water, before assembly, disassembly or cleaning, always make sure to unplug the appliance from the el.

- For increased safety, we recommend that you install a residual-current device in the bathroom's el. circuit (RCD) with rated residual operating current not exceeding 30 mA. Consult with an inspecting technician or electrical engineer.

NOTICE

- The appliance must only be powered by a safe low voltage of 5 V.
- For power supply, use only a USB source complying with valid standards with an output voltage of 5.0 V and an output current of at least 1000 mA (a mobile phone charger with a double-insulated USB connector or a backup power source is suitable).
- The appliance does not have its own power supply.
- Do not use USB ports on other devices (e.g. computer, laptop, tablet, mobile phone, etc.) for power supply. These devices may not be sized for the current consumption of the appliance and could be damaged.

SPECIAL SAFETY INSTRUCTIONS

- The appliance should not be used by persons with reduced sensitivity of the oral cavity. It may lead to injuries or damage teeth and gums. If you have any doubt about possible use in your specific case, please consult your dentist before using the appliance.
- The appliance is not to be used by children, incapacitated persons or persons unable to express their sensation.
- If you have undergone tooth or gum surgery recently, suffered from long-lasting tooth disorder or have undergone extensive treatment in the last two months, or if you suspect any tooth-related problems (fillings, crowns, etc.), please consult your dentist before using the appliance.
- Should any health problems occur (e.g. extensive bleeding, or if bleeding does not stop after 2 weeks of use), stop using the appliance and ask your dentists for advice.

GENERAL SAFETY INSTRUCTIONS

- The appliance is a household appliance intended for individual use for personal dental care! It is not designed for commercial use, or to be used by several patients in dental offices or other medical facilities!
- Do not plug or unplug the power cord to/from the electrical socket with wet hands!
- Never immerse the adapter into water or other liquids (not even partly) and do not wash it in running water!
- The appliance must never be left running unsupervised!
- Keep the appliance in a dry place.
- Use only accessories intended for this type of appliance. Do not use damaged / deformed attachments.
- Do not use the application in exterior or where aerosol products are used (such as sprays), or where oxygen is handled.
- For hygienic reasons, one nozzle must not be shared by more people. Otherwise, infections or inflammation in the oral cavity might be transmitted.

- Do not use/recharge during bathing.
- Keep the accessories (i.e. nozzles) out of reach of children and incapacitated persons.
The nozzle is small enough to get stuck in your throat or be swallowed. If anyone swallows a nozzle, they need to seek immediate medical attention!
- Do not press the nozzle tip against teeth or gums with excessive force.
- Do not put the appliance in operation if the reservoir is empty or if the nozzle is removed.
- Do not insert the appliance or its accessories in any body cavity other than the oral cavity.
- Do not direct the water jet under the tongue, in the throat, ear, nose, eyes or other sensitive areas as it might result in injuries.
- Never direct the water jet at people, animals, plants or appliances containing electrical parts, or at their power cords.
- The appliance is not designed for animal dental care.
- Only use the appliance in a horizontal position in places where it cannot overturn, at a sufficient distance from sources of heat (e.g. stoves, fireplaces, or thermal radiators), and sources of increased air humidity (e.g. sauna, pool).
- Protect the appliance from dust, chemicals, excessive moisture, toiletries, liquid cosmetics, large temperature changes and do not expose it to direct sunlight.
- As a result of production testing, dew might form on the inside of the suction hose or reservoir. Do not consider this to be a fault when purchasing the product.
- Do not overfill the reservoir. The maximum quantity of water is 250 ml.
- Never fill the reservoir with hot water (max. 40 °C), mouthwash, essential oils, water enriched with alcohol, salt, vinegar, descaling agents, softeners or other cleaning chemicals. Never put any items in the reservoir.
- The appliance must be stored at a temperature over 0 °C and it needs to acclimatize to a room temperature before use.
- Always clean the appliance after use. Failure to do so might result in deterioration of your state of health as fungi and bacteria might grow unless you disinfect the appliance.
- Never try to repair the appliance yourselves.
- Do not let the power cord hang loosely over the edge of the worktop where children might get hold of it.
- Check the condition of the power cord regularly.
- The appliance gets hot when charging, which is absolutely normal.
- Charge the battery at a temperature close to the room temperature.
- Do not dispose the battery into fire. There is a risk of explosion!
- Do not charge a battery leaking electrolyte.
- Electrolyte leak is caused by overloading the appliance or using it at extremely high temperatures. If you spill electrolyte on yourself, wash the area with water and soap and rinse with lemon juice with vinegar. If electrolyte gets in your eyes, flush with clear water for several minutes and call a physician immediately.
- Do not expose the battery to temperatures higher than 50 °C to prevent battery damage.
- Never charge the appliance in direct sunlight or close to heat sources, and if the temperature is below 10 °C or above 40 °C.
- The power cord must not be damaged by sharp or hot objects or open flame, it must not be submerged in water or bent over sharp edges.
- Switch the appliance off immediately after use and put it in a dry place, out of reach of children and incapacitated persons.
- Never use the appliance for any purpose other than the intended purpose described herein!
- The manufacturer shall not be held liable for any damage caused by improper handling of the appliance and its accessories (e.g. injury, damage to the appliance, fire, etc.)

II. DESCRIPTION OF THE APPLIANCE AND ACCESSORIES (pic. 1)

A – Oral irrigator body

- A1 – ON/OFF switch
- A2 – MODE switch / LOCK button
- A3 – active mode indicator C, S, P, D
- A4 – nozzle ejector button
- A5 – charging socket

- A6 – charging socket cover
- A7 – water reservoir
- A8 – reservoir plug
- A9 – lid lock
- A10 – lid

B – Nozzle for teeth and gums (2 pcs included)

C – Orthodontic nozzle for fixed braces

D – Power cable

III. INSTRUCTIONS FOR USE

Remove all packaging materials, take out the oral irrigator and any accessories. Put the appliance on a flat and dry surface (so that it cannot fall into water), out of reach of children.

CHARGING

After removing the cover **A6** plug the power cord into the **A5** charging socket. The other end to a source such as a compatible AC power adapter or backup power source.

A discharged battery is signaled by the flashing of the **A3** indicators 5 times. Charging is indicated by gradual flashing of the **A3** lights; full charge then the steady light of these lights. After charging, replace cover **A6**. For a full charge, the battery needs to be charged for approximately 3 hours. Depending on the cleaning mode, can empty approximately 40 – 45 water reservoirs **A7** (about 50 min.), which will last for about 2 – 3 weeks if used twice a day for two minutes as well.

Note: Although the brush / oral irrigator should not be charged before they are fully discharged, we recommend charging them at regular intervals in order to ensure maximum battery life (e.g. after 2 weeks of use).

Caution

- Do not exceed the recommended time of charging 3 hours during regular charging!
 - The appliance cannot be used while charging.
 - Use only the type of battery approved by the manufacturer in the appliance.
- Do not use batteries that must not be recharged.

IV. USE

The water system must be flushed out before first use. Fill the water reservoir **A7** with clean drinking water according to the following instructions, re-assemble the appliance and switch it on. Allow all the water to run through the oral irrigator to the washbasin. Repeat this procedure at least twice.

Filling with water

It is necessary to open the container plug **A8** first and reinsert it properly after filling. See Figure 2. The bin **A7** can be opened by pressing the lock **A9** and pulling out the lid **A10**.

Attachment B/C – the snozzle into the opening at the top of the oral irrigator body A until it locks in place (you will hear a click). “Pull” the attachment slightly so as to check that it is properly fitted. To remove the attachment, press the ejector button A4 first and take out the attachment pulling it gently.

Button control

Using the mode switch button **A2**, select a cleaning mode first:

- a) **Clean mode** – standard intensive mode for regular cleaning of healthy teeth and gums,
- b) **Soft mode** – gentle mode for sensitive teeth and gums,
- c) **Pulse mode** – pulse mode for gentle gum massaging.
- 4) **mode D (DIY)** – allows you to adjust the intensity yourself. After activating this mode and turning on the appliance, press and hold button **A2** - the intensity of the water jet will be switched automatically. You can set the desired intensity by releasing the button at the moment when the intensity that suits you is active.

The selected mode can be changed even when cleaning is in process. Activation of the mode selected is indicated by the indicator **A3** coming on. Activate the appliance by pressing the ON switch **A1**. When you press the ON switch again, the appliance will switch off or after the 2 minutes expires, the appliance will automatically turn off. After switching off/exiting the mode, the number of indicator lights **A3** corresponding to the current battery charge status (4 = max. charge, 1 = min. charge) will remain lit for a few moments.

Note

- The appliance remembers the last set mode.
- After emptying the container, switch off the appliance and, if necessary, refill it.

Method of use (Fig. 3)

As plaque which clings to the teeth surface cannot be removed simply with a stream of water, it is necessary to clean your teeth thoroughly with a toothbrush before using the oral irrigator. As soon as you have cleaned your teeth replace the brush with the oral irrigator. Attach the nozzle **B** to the oral irrigator and put the nozzle tip between two teeth at the gum line level at a 90 degree angle. To prevent the water from splashing, press the lips slightly together around the nozzle, but keep the mouth slightly open so that water can flow out smoothly. Lean over the washbasin while cleaning your mouth. Slowly guide the nozzle until you feel it settle between the next two teeth, smoothly follow the gum line and the shape of your teeth. Continue doing so until you clean all inter-dental spaces, including the molars, in all areas of the mouth. Do not apply any pressure on the nozzle, let it do all the work on its own. Depending on the selected mode of cleaning, the full water reservoir **A7** will last for about 60 – 90 seconds. For optimal mouth irrigation, we recommend refilling the reservoir 2 times.

When cleaning is finished

- Open the water reservoir **A7** and empty it, pouring any remaining water in the washbasin.
- Switch the oral irrigator on for a few seconds so as to empty the inside as well.
- Wipe all the components with a dry cloth.

Buttons lock

Long press the **A2** button in the off state to activate the button lock against accidental activation (e.g. when travelling). Activation is signaled by the lighting of all **A3** lights. You can deactivate the lock in the same way.

Caution

- To avoid proliferation of fungi and bacteria if the appliance is not used for a longer time, empty the reservoir thoroughly after each use!
- During the first few days of using the oral irrigator, your gums might bleed slightly but this usually stops on its own after a few days.
- To ensure the best results when cleaning, we recommend replace every 6 replace every 6 months the nozzle. However, if the attachment are showing signs of significant wear or damage, replace them sooner.

V. MAINTENANCE

For maximum hygiene, clean the oral irrigator after each use! Failure to do so might result in deterioration of your state of health due to proliferation of fungi and bacteria. Do not use any rough and aggressive cleaning agents! The nozzle for teeth and gums **B**, **C** under running water. Use the damp cloth to wipe the oral irrigator body **A**, cable **D** as well. For hygienic purposes we recommend that you dip the nozzle / tongue scraper into a denture cleaner every 2 – 4 weeks.

Alternatively, you can use other appropriate disinfectants. When using cleaning agents / disinfectants, always follow the instructions indicated on the packaging.

Make sure that the contact surfaces and sealing elements are clean, not damaged and properly positioned. If you are not going to use the appliance for an extended period of time, empty the entire water system thoroughly (i.e. empty the reservoir, press the switch A1, drain the water until it stops flowing from the nozzle).

Other important notes:

- Do not clean the appliance in a dishwasher and do not use any aggressive cleaning agents for cleaning!
- Never use water at a temperature exceeding 50 °C to clean the oral irrigator
- Allow all the appliance components to dry properly before their re-use.

Store the appliance and all the accessories in a dry, dust-free and safe place, out of reach of children and incapacitated persons.


Replacement attachment **B** can be purchased under order numbers ETA470890100.

Replacement attachment **C** can be purchased under order numbers ETA470890200.

Specialized service is required for overhaul maintenance or maintenance requiring disassembly of the internal components of the appliance! The right for warranty repair expires if the above mentioned instructions are not adhered to.

More extensive maintenance or maintenance requiring intervention in the inner parts of the appliance must be carried out by professional service! Failure to comply with the instructions of the manufacturer will lead to expiration of the right to repair!

The information regarding warranty and product repairs you can find on the website www.eta.cz/servis-eu.

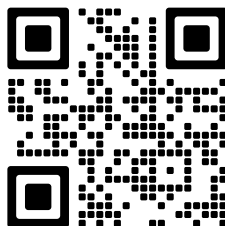VI. ENVIRONMENT

If the dimensions allow, there are labels of materials used for production of packing, components and accessories as well as their recycling on all parts. The symbols specified on the product or in the accompanying documentation mean that the used electrical or electronic products must not be disposed of together with municipal waste. For proper disposal, hand them over at special collection places where they will be accepted free of charge. Suitable disposal of the product can help to maintain valuable natural resources and to prevent possible negative impacts on the environment and human health, which could be the possible consequences of improper waste disposal. Ask for more details at the local authorities or at a collection site. Fines can be imposed for improper disposal of this type of waste in agreement with the national regulations. Dispose of the battery safely by means of collection networks intended for this purpose. Never dispose of the batteries by burning!

Replacing battery

Remove the batterie from the appliance only when they are fully discharged. Screw out all bolt and open the pedicure set. Gradually disconnect and remove the batterie. Proceed as in the picture 4.

VII. TECHNICAL DATA

Voltage (V)	shown on the type label of the appliance
Weight (kg) approximately	0.25
Appliance Protection Class	III.
Size, (mm):	65 x 75 x 210
Operating pressure (mode C,S,PD)	30-120 PSI

The manufacturer reserves the right to change the technical specifications and accessories for the respective models.

The capacity of the product battery decreases over time. Decreasing of battery capacity occurs during normal use of the appliance and is considered usual wear and tear of the appliance. Decreasing battery capacity is not a defect of the appliance.

WARNINGS AND SYMBOLS USED ON THE APPLIANCE, PACKAGING OR IN THE INSTRUCTIONS MANUAL:


HOUSEHOLD USE ONLY. DO NOT IMMERSE IN WATER OR OTHER LIQUIDS.

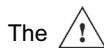TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, CARRIAGES OR PLAYPENS. THIS BAG IS NOT A TOY.

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© DATE 05/05/2023
Service: Omega electric GmbH, Servicezentrale, Industriering 2, D- 04626 Schmölln. Bei Fragen zu