




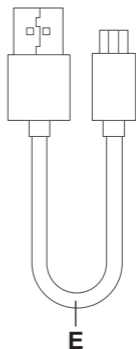
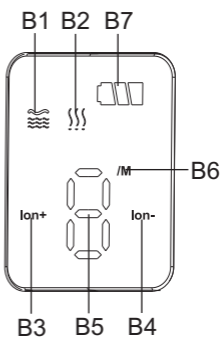
<i>Kosmetická ultrazvuková špachtle s galvanizací • NÁVOD K OBSLUZE</i>	<b>CZ</b>	<b>5-10</b>
<i>Kozmetická ultrazvuková špachtľa s galvanizáciou • NÁVOD NA OBSLUHU</i>	<b>SK</b>	<b>11-16</b>
<i>Ultrasonic spatula with galvanisation • USER MANUAL</i>	<b>EN</b>	<b>17-22</b>
<i>Kozmetikai ultrahangos spatula galvanizálással • KEZELÉSI ÚTMUTATÓ</i>	<b>HU</b>	<b>23-28</b>
<i>Ultradźwiękowa szpachelka kosmetyczna z galvanizacją • INSTRUKCJA OBSŁUGI</i>	<b>PL</b>	<b>29-34</b>
<i>Ultraschallspachtel mit der Galvanisierung • BEDIENUNGSANLEITUNG</i>	<b>DE</b>	<b>35-41</b>
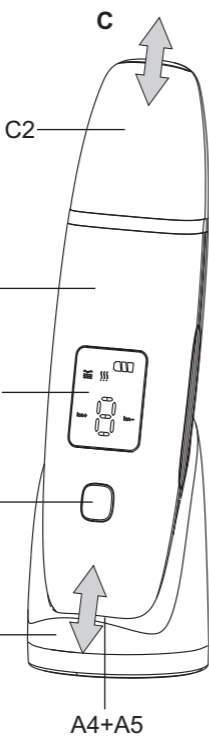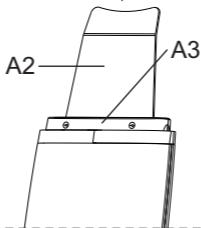
# Fenité

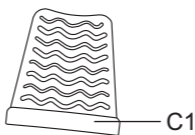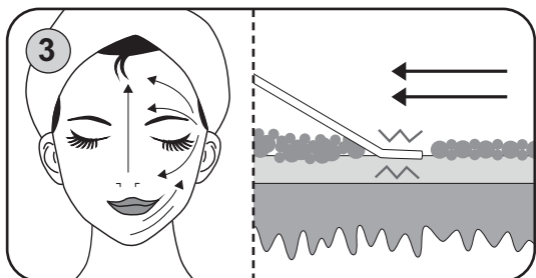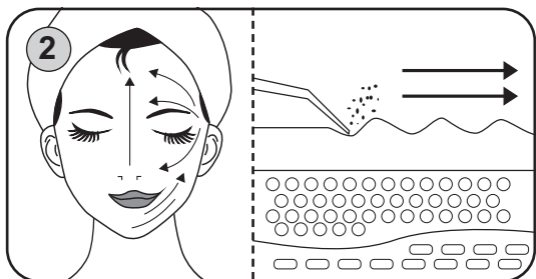1

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*TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, CARRIAGES OR PLAYPENS. THIS BAG IS NOT A TOY. Nebezpečenstvo udusenía. Nepoužívajte*

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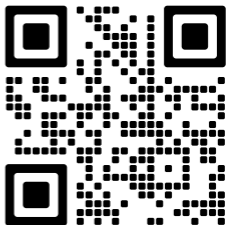
Fenité

## INSTRUCTIONS FOR USE

Dear customer, thank you for purchasing our product. Please keep these instructions including the receipt and, if possible, the carton box with the internal packing.

**I. SAFETY WARNING**

- Please consider the instructions herein as a part of the appliance and pass them over to the next user of the appliance.
- The appliance is only designated to be used in household and for similar purposes. It is not designed to be used in beauty parlours or for any other commercial utilisation !
- **This appliance may only be used by children aged 8 years or older, and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.**
- **Never immerse the appliance in water!**
- **When handling the appliance, please proceed so that no personal injury or damage to the appliance will occur.**
- **Never use the appliance if it does not correctly operate, fell on the ground or suffered any damage.**
- Protect the appliance from dust, chemicals, excessive humidity, toiletries, aggressive liquid cosmetic preparations and intensive temperature changes, and do not expose it to direct sunlight.
- Before use, take off all metal items and jewellery (e.g., glasses, watch, earrings, bracelets, rings, etc.).
- Use only the type of battery approved by the manufacturer in the appliance. Do not use batteries that must not be recharged.
- Do not insert the appliance or its accessories in any bodily orifices.
- Never use the appliance for other than the intended purposes as described in these instructions for use!
- Only use the appliance with the accessories designed for this type. Use of unapproved accessories may cause injury.
- **WARNING:** In the case of incorrect use of the appliance that contradicts the operating instructions, there is a risk of injury.
- Texts in a foreign language and images provided on the packaging or on the product, if any, are translated and explained at the end of this language version.
- This appliance contains batteries that are not replaceable.

- The manufacturer is not responsible for any loss caused by improper handling of the appliance and its accessories (e.g. injury, harm to health) and is not responsible according to the appliance in the case of non-compliance with the safety instructions above. This appliance is not intended as a substitute for medical care.
- Do not use the appliance if your skin is oversensitive, damaged or exematic (e.g. in the case of fresh scars, small injuries, burns, erythema, swollen cheeks, skin diseases...).
- If you feel pain or experience any unpleasant feeling during the cleaning, or if any other health problem appears, immediately switch off the appliance, and consult its further use with your practitioner or dermatologist.
- Do not use the appliance with exfoliating or any other aggressive gels (for example with alcohol)!
- The appliance is not intended to cure any skin disease symptoms. Please consult use of the appliance with your doctor, where necessary.
- Do not use the appliance on children, incapacitated persons or persons unable to express their feelings.
- Do not use the appliance near eyes or other vulnerable body parts.

#### SPECIAL SAFETY INSTRUCTIONS

- Children must be kept away from the appliance and its adapter.
- The appliance has no charging source of its own, it is only designed for charging via USB (an adapter with USB output, PC with USB output with output voltage of 5.0 V, and output current of at least 1000 mA). Before making any connection, please pay attention to instructions of the manufacturer of the respective power source.
- Never use the appliance if the adapter is damaged, if the appliance does not operate properly, fell on the ground and suffered any damage, or fell into water. In such case bring the appliance to the specialised service centre to verify the appliance safety and proper operation.
- For additional protection we recommend installing a protective switch (RCD) in the electrical circuitry of the bathroom, with the nominal actuating current not exceeding 30 mA. Ask an inspecting technician or an electrician for advice.
- The appliance and the connected adapter must not be used at places

where it could fall into bath tub, wash basin or swimming pool. If the appliance fell into water anyway, do not take it out! First of all, unplug the adapter from the socket and only then take the appliance out. Then bring the appliance to a specialised service centre to check its safety and proper functionality.

- To power the appliance, only apply safe small voltage in accordance with the type plate.
- If the appliance is used in a bathroom, it must be disconnected from the mains after use by unplugging the power cord from the wall socket, because the appliance poses a risk if near to water, even when switched off.
- Always unplug the appliance from charger if you leave it unattended and cleaning.
- The mains adapter supplied with the appliance is only designated for this appliance.
- **WARNING:** The appliance is not suitable for use when wet (e.g. in the shower, etc.).
- Never plug the adapter in the electrical socket or unplug it with wet hands.
- Never immerse the adapter into water or other liquids (not even partly) and do not wash it in running water! At the same time, prevent them from being splashed with liquids.
- The appliance gets hot when charging, which is absolutely normal.
- Charge the battery at a temperature close to the room temperature.
- Never charge the appliance in direct sunlight or close to heat sources, and if the temperature is below 10 °C or above 40 °C, adversely affects the accumulator service life.
- Prevent hanging the power adapter cord loosely over desktop edge where children may reach it.
- The power cord of the adapter must not be damaged with sharp or hot objects or open fire, and must not be immersed into water or bent over sharp edges.

## II. DESCRIPTION OF THE APPLIANCE (Fig. 1)

### A – Appliance body

- A1 – ON/OFF button + operating mode setting button
- A2 – Stainless steel applicator
- A3 – Blue light photon therapy
- A4 – USB charging socket
- A5 – Contacts for the charging stand

**B – Backlit LCD**

- B1 – Set mode icon - high intensity
- B2 – Set mode icon - low intensity
- B3 – Set mode icon - positive ionisation
- B4 – Set mode icon - negative ionisation
- B5 – Remaining operation time (in minutes)
- B6 – Symbol Minutes
- B7 – Battery level indicator

**C – Accessories**

- C1 – Silicone case (skin cleansing)
- C2 – Protection case

**D – Charging stand****E – USB cable \***

\* USB cable: A source with USB terminal may be connected to USB cable (for example AC power adapter, backup power source, etc.).

## III. INSTRUCTIONS FOR USE

Before using the appliance, its battery must be charged first. Charging is possible with the enclosed USB cable from the mains, i.e. an AC adapter, or a backup power source. Charging stand **D** can also be used for charging the battery.

**CHARGING**

- Before the first use or if the appliance was not used for a long time, charge it completely.
  - The standard charging time is **approximately 2 hours**, and the fully charged accumulator allows for the operating time of **up to 100 minutes**. The operating time decreases depending on the appliance loading and due to the accumulator ageing and wear and tear.
1. Make sure the appliance is switched off.
  2. Connect the appliance to the power source (battery charger) using USB cable or use the charging stand, which can also be connected to the power source (battery charger) using USB cable, and then place the appliance on the stand correctly (as shown on fig. 1). If the appliance is turned on and you start charging, the appliance shall switch off.
  3. Do not exceed the designated charging time!

**Note**

Battery level and charging state are displayed by indicator **B7** (the indicator is flashing during the charging process and lighting continuously as soon as the charging process has been completed).

**Optimise the battery service life**

- Do not leave USB cable connected to the appliance.
- During charging or regular operation, the appliance heats up, which is a normal phenomenon, and not a reason for a complaint.

**CONTROL**

- By long pushing of button **A1** (approx. for 1 sec.) you switch the appliance on/off. Upon switching on and off an acoustic signal sounds.
- Press short **A1** button to switch over the modes. Currently set mode is shown by LED indicators on the display (**B1, B2, B3, B4**).

**USE**

With use of acoustic vibrations (frequency of 24 000 Hz) the appliance thoroughly cleans the skin and removes dead skin cells. It also improves absorption of various cosmetic preparations and softens the skin.

Galvanic current technology:

- positive ionisation (Ion+) – for deep skin cleansing
- negative ionisation (Ion-) – for deep application of cosmetics

The appliance also includes photon therapy function using blue light (antibacterial and anti-inflammatory effects).

#### **Skin cleansing (Fig. 2)**

1. Moisten the facial skin.
2. Switch on the appliance and set the mode of operation you want.
3. Apply **A2** stainless steel applicator on the face. Lightly pressing, move the applicator evenly on the entire facial skin, as shown on Fig. 2.
4. Switch off the appliance after application and rinse the skin with water.

#### **Note**

Use silicone case **C1** for fine cleansing and skin massage.

#### **Application of cosmetic preparations (Fig. 3)**

1. Apply the cosmetic preparation on the face.
2. Switch on the appliance and set the mode of operation you want.
3. Apply **A2** stainless steel applicator on the face. Pressing the applicator lightly on the skin, move it evenly on the entire facial skin, as shown on Fig. 3., allowing the preparation to be absorbed through the skin.
4. Switch off the appliance after application.

#### **Auto-off mode**

After approximately 9 minutes of continuous operation the appliance automatically switches off. Before switching off, an acoustic signal sounds. The remaining operation time is displayed by indicator **B5**.

## **IV. MAINTENANCE**

The appliance must be cleaned after every use! Use only soft and damp cloth to wipe the surface of the appliance, charging stand and USB cable. Stainless steel applicator **A2** can also be wiped by soft and damp cloth. Silicon case **C1** for skin cleansing can be rinsed in flowing water. Allow to dry properly before using it again. Do not use any abrasive or aggressive cleaning agents!

**The capacity of the product battery decreases over time. Decreasing of battery capacity occurs during normal use of the appliance and is considered usual wear and tear of the appliance. Decreasing battery capacity is not a defect of the appliance.**

## **V. ENVIRONMENT**

If the dimensions allow, there are labels of materials used for production of packing, components and accessories as well as their recycling on all parts. The symbols specified on the product or in the accompanying documentation mean that the used electrical or electronic products must not be disposed of together with municipal waste. For proper disposal, hand them over at special collection places where they will be accepted free of charge. Suitable disposal of the product can help to maintain valuable natural resources and to prevent possible negative impacts on the environment and human health, which could be the possible consequences of improper waste disposal. Ask for more details at the local authorities or at a collection site. Fines can be imposed for improper disposal of this type of waste in agreement with the national regulations. Dispose of the battery safely by means of collection networks intended for this purpose. Never dispose of the batteries by burning!

#### **Taking out the battery**

Only take out the battery from the appliance after it is fully discharged (i.e. let the appliance in operation until its motor stops), and after the end of the appliance's service life! The battery should not be removed and is not accessible during normal operation! With use of a suitable tool demount the appliance body and replace the accumulator. However pay attention to prevent any cut injury.

## VI. TECHNICAL DATA

Voltage (V)	Specified on the rating place
Appliance weight (g)	94
Appliance protection class	III.
Dimensions (H x D x W) in mm	20 x 168 x 49

**The manufacturer reserves the right to change the technical specifications and accessories for the respective models.**

*HOUSEHOLD USE ONLY. DO NOT IMMERSE IN WATER OR OTHER LIQUIDS.*

*TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, CARRIAGES OR PLAYPENS. THIS BAG IS NOT A TOY.*

The  symbol indicates a WARNING.

Do not use the appliance close to bath tubs, showers, wash basins or other water containers.


Washable under running water

## VII. CUSTOMER CARE AND SERVICE

**More extensive maintenance or maintenance requiring intervention in the inner parts of the appliance must be carried out by professional service! Failure to comply with the instructions of the manufacturer will lead to expiration of the right to repair!**

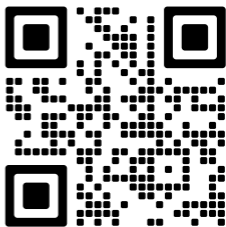The information regarding warranty and product repairs you can find on the website [www.eta.cz/servis-eu](http://www.eta.cz/servis-eu).

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© DATE 09/06/2022