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Analytic bathroom scale with a Smart application • INSTRUCTIONS FOR USE **EN** 23-30

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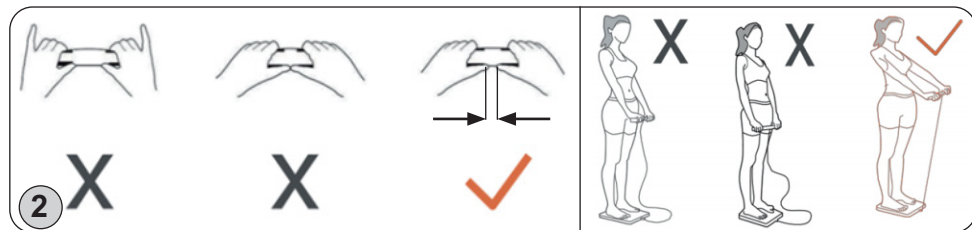
EN: when measuring, it is necessary to use the correct standing posture for measurement: pull up the handle, keep your back and knees straight, lift your hands horizontally to your chest, and straighten your elbows. both hands and the handle electrodes need to be completely sealed, and both hands must not interfere with each other.

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TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, ARRIAGES OR PLAYPENS. THIS BAG IS NOT A TOY.

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Analytic bathroom scale with a Smart application

eta 8779

Vital Condition

INSTRUCTIONS FOR USE

Dear customer, thank you for purchasing our product. Please read the operating instructions carefully before putting the appliance into operation and keep these instructions including the receipt and, if possible, the box with the internal packing.

Tip

- For initial pairing with your home Wi-Fi network using the “ETA Body Balance” app, please visit our website www.etasince1943.com/care.

I. SAFETY WARNING

- Consider the instructions for use as a part of the appliance and pass them on to any other user of the appliance.
- **The appliance is intended for use in households only and for similar use! It is not intended for use in health-care facilities and for commercial use!**
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- Never use the appliance if the USB power cable or adapter is damaged, if it does not work properly, if it fell to the ground and got damaged. In this case take the appliance to a special service to check its safety and proper function.
- The appliance must only be powered by a safe low voltage of 5 V.
- The appliance has no charging source of its own, it is only designed for charging via USB (an adapter with USB output, with output voltage of 5.0 V, and output current of at least 0,5 A). Before making any connection, please pay attention to instructions of the manufacturer of the respective power source.

WARNING

- **To power the appliance, only use a source of safe low voltage of 5 V which complies with applicable standards (mobile phone chargers with a USB-C connector with double insulation are suitable).**
- **Do not use USB ports on other devices (e.g. computers, laptops, tablets, mobile phones, etc.) to power the appliance. Such devices may not be dimensioned for the current consumption of the appliance and it.**

- Never use the appliance if it malfunctions, if it fell to the ground and got damaged or if it fell into water. In these cases bring the appliance to a specialised service to have the safety and functions of the device checked.
- **We recommend: Do not step onto the scale with wet feet or on wet surface of the scale, there is a risk of slipping.**
- **Never immerse the scale into water (even their parts)!**
- Keep in mind that the data read from this scale are approximate and inaccurate in comparison to results of actual medicinal analyses. The exact determination of biometrical characteristics can be done only by a doctor via medical methods.
- Do not use during pregnancy.
- Do not stand on the edge of the scales!
- Do not use the fat analysis feature on this scale if you are using one of the following devices: pace maker, electrocardiograph and/or other medical device implanted in the body or device used as support.
- If the scale was stored in lower temperatures, let it adjust to the room temperature.
- Use the scale in a horizontal position only and at places where they cannot knock over, in sufficient distance from heat sources (e.g. **stove, fireplace, heat radiators**), sources of increased air humidity (e.g. **sauna, bathroom, swimming pool**) and equipment with strong electromagnetic field (e.g. **microwave oven, radio**).
- Protect the personal scale from dust, chemicals, excessive humidity, toiletries, liquid cosmetic preparations, intensive temperature changes and do not expose them to direct sunlight.
- Weigh yourself without clothing and footwear, before meals and always at the same time of day.
- Do not place the scale on uneven surface or on surface covered with carpet. Unstable or soft surface under the scale can have negative effect on precision of weighing.
- In handling the scale, make sure injury of persons or damage of the appliance is prevented.
- Do not overload the appliance weight greater than the maximum weight capacity! This prevents damage to the appliance.
- Handle the scale and the handheld sensor and its connecting cable with care. Do not throw it or jump on it. Do not pull the hand-held sensor with force to avoid damaging the connecting cable and winding mechanism. The scale is durable, but improper handling can destroy its electronic sensors (fig. 2).
- Do not put any objects on the surface of the scale during storing, they could be damaged.
- To power the appliance, only apply safe small voltage in accordance with the type plate.
- Do not leave the appliance plugged in power supply without supervision!
- Always unplug the appliance from charger if you leave it unattended or before cleaning.
- **Do not put or take out the plug to or from a socket with wet hands!**
- Charge the accumulator at the normal room temperature. Do not expose the accumulator to temperatures higher than 40 °C to prevent causing damage to the accumulator.
- The appliance must be disconnected from the power adapter when removing the battery.
- Do not through the battery into fire. There is a risk of explosion!
- This appliance contains batteries that are not replaceable. Only remove the accumulator if the appliance is to be discarded (see article **IX. ECOLOGY**).
- The power cord of the adapter must not be damaged with sharp or hot objects or open fire, and must not be immersed into water or bent over sharp edges.
- Regularly check the condition of the charging cord.
- Never use the appliance for any other purpose than for that described in these instructions for use!
- **WARNING:** There is a risk of injury in case of incorrect use (not in accordance with the manual) of the appliance.

- The producer is not responsible for damage caused by improper use of the appliance (e.g. **breaking the board, etc.**) and it does not cover the appliance in the case of failure to comply with the safety warnings above.


Important notes

- This scale requires you to weigh yourself barefoot for the reason of correct conductivity during the body fat analysis.
- Do not place the scale on an uneven surface or a surface covered with a carpet.
- Excessive drinking, eating, exercising, health conditions, used medicaments, women's menstruation cycle etc. influence the results of the analysis.
- There must not be a contact between feet, calves or thighs. Otherwise the measurement could be inaccurate.
- Due to the inherent differences between people, the instructions and the application are for information purposes only.
- The hand-held sensor must also be used for the analysis to work properly. Simply stepping on the weighing surface is not enough. Grasp the sensor with the whole palm of your hand as shown in the figure below. Grasp the left arm with the left palm and the right arm with the right palm.

Wireless connectivity and compatibility

- ETA can in no way be held liable for any accidental, indirect or subsequent damage, data loss or any other loss incurred as a result of information leakage during communication through BLUETOOTH, Wifi technology.
- 100% successful wireless communication cannot be guaranteed with all Smart devices. Devices have to be compatible with BLUETOOTH specifications set forth by the company Bluetooth SIG, Inc., and their compliance with the specifications has to be certified. However, given the variety of devices available in the market, even if a device is in compliance with the specifications, there might be cases when the character or specifications of the given Smart device do not allow connectivity, or hinder connectivity in any way. This is not a reason to make a warranty claim.
- Radio waves transmitted by devices may interfere with operation of medical equipment. Such interference may cause a failure. Therefore, do not use the device in the following places:
 - in hospitals, trains, planes, petrol stations and in areas where flammable gases may be present.
 - in the vicinity of automatic doors or fire detectors.

The application is owned and operated by a third party (hereinafter the "Subject") different from HP TRONIC Zlín, Ltd. (hereinafter the „Seller“). No personal information is transmitted between the Subject and the Seller. When registering for the application, the customer transfers the personal data directly to the Subject, which becomes the personal data controller. The Application and the associated privacy policy are governed by the Rules of the Subject. Development and any changes to the application are made solely by the owner of the application.

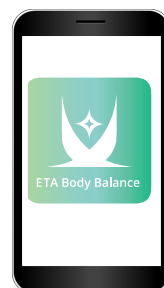II. GENERAL INFORMATION

You have purchased an advanced device that is a multifunctional modern assistant for monitoring and analyzing various body characteristics. The handheld sensor makes the results even more accurate. Together with the **ETA Body Balance** app (compatible with Android or iOS), which you connect to the scale via Bluetooth, you will have the results clearly in one place and available from any compatible device. After the scale is paired with the Smart device, the **ETA Body Balance** app starts displaying the measurement results in real time. The scale itself only displays the weight, but the application also receives other measured values for further processing. Before first use, download the **ETA Body Balance** application to your Smart device for free. You can find the application in the Apple Store or Google Play database, or to speed up the process, you can use the following QR codes:

Note

- If the QR code is not active, find the application manually by name.
- All other information, operating instructions of the app, instructions and hints for successful pairing with home network can be found at: www.etasince1943.com/care.

The application requires an active internet connection.

III. FEATURES OF THE DEVICE (FIG. 1)

A – Unit of scale

- | | |
|--|-----------------------|
| A1 – handheld sensor | A4 – power USB-C slot |
| A2 – display | A5 – cover |
| A3 – buttons – (reset to factory settings) | |

B – USB-C cable (a source with USB terminal may be connected to USB-C cable (for example AC power adapter, backup power source, etc.)

IV. PREPARATION FOR USE

Remove all the packing material and take the scale out. Remove all possible adhesion foils, stick-on labels or paper from the scale. Place the scale on flat, stable and dry surface.

Charging of accumulator (Fig. 3)

Before the first use or if the appliance was not used for a long time, charge it completely. Connect the corresponding USB-C cable to the power supply terminal **A4**, and connect the other side to a power source, for example to a compatible AC power adapter, or to a backup power source (please confirm the voltage of the external electric device is =5 V / 0,5 A) before charging it. The first charging or charging after a long period of time is approx. 120 minutes! The operation of charging of the accumulator is indicated by a display **A2**. Unplug the adaptor from the wall socket first and only then unplug the connector from the slot. Do not exceed the charging time given.

ATTENTION

- Do not leave the USB-C cable permanently plugged in the appliance.
- When the battery of the accumulator is low and it needs to be recharged, „**Please charge in time**“ appears on the display.
- If you're charging the accumulator for the first time after being stored away for some time, the accumulator will charge to about 60 % of total capacity. After discharging and recharging again, the accumulator will charge to 100 %.
- To extend the lifespan of the accumulator, we recommend discharging the accumulator completely 2x a year by normal use and then charging it fully again.

Weight and height units

Weight units can also be changed in the application.

kg = kilogram	1 kg = 2.2 lb
lb = pound (British weight unit)	1 lb = 0.454 kg
st = stone (American weight unit)	1 st = 6,356 kg
Inch = inches (British length unit)	1 In = 2,54 cm

V. INSTRUCTIONS FOR USE

Your scale is fine electronic equipment. One of the things that are most difficult to weigh is human body as it is always in motion. If you want to obtain reliable results, always try to step with your feet at the same place on the scale and stand still. If possible, weigh yourself at the same time of a day (preferably in the morning) to obtain comparable results. Remember that excessive clothes can increase your weight as well as eating and drinking before weighing. Short-term weight fluctuation is generally attributed to loss of liquids. The handheld sensor, which is a fixed part of the scale, adds an additional element to ensure more accurate analysis results.

Scale initialization (calibration), (Fig. 3)

This process needs to be repeated whenever you move the scale.

- Step up on the centre of the weighing scale which turns it on and then step down.
- “0.0” will appear.
- Wait till it switches off automatically. The scale is now ready for an analysis.

1) STANDARD WEIGHT MEASUREMENT

Without the use of the **ETA Body Balance** app in the Smart device, the scale can provide basic weight measurements as any other scale.

Gently step on the scale (necessary with a handheld sensor) - make sure your palms are fully gripping the electrodes, then the scale will be automatically turned on.

Stand evenly on the scale without moving and wait until your weight shown on display.
 – After a moment, the data will be automatically displayed on the screen. **Then the scale are switched off automatically.**

Note

- The scale is switched on automatically if you do not step on it within approx. 10 seconds after „0.0“ is displayed.
- We recommend using the scale from about 10 kg of the subject.
- The icon flashes on the display to indicate communication / searching for nearby devices. When the icon stays lit, the connection with the application is on.
- Keep the app turned off for normal weighing.

2) WEIGHT MEASUREMENT WITH AN ANALYSIS AND USE OF THE APP

Download the app, create an account and pair the app with the scale. Turn the app on and then weigh yourself barefoot in the standard way (necessary with a handheld sensor) - make sure your palms are fully gripping the electrodes. After a successful measurement, the data is automatically transferred to the app.

Pairing the scale

Follow this procedure both for the first pairing and if you need to pair during regular use. Turn on the Bluetooth function on your Smart device.

Turn on the ETA Body Balance app,

The first time you turn it on, the basic working interface will be displayed. Click on „Add Device“ in the top right corner.

Turn on the scale.

Once the Smart device has found the scale, it will show it on the display.

For more information, download the instructions for use of the application, see chap.II.

VI. TROUBLESHOOTING OF POSSIBLE PROBLEMS

If you have problems with weight:

- Check whether the scale is on a flat floor and they do not touch the wall.
- If nothing appears on the scale when you want to switch them on by tapping on them or if „**Please charge in time**“ appears only, appears only, charge the battery.
- If „**Overload**“ appears on the display, the scale was overloaded.
- If signs appear, such as if the display is flashing or if there is a faulty sign, try to move the scale away from the influence of a source of interference or ensure switching off the source for the time of using the scale.
- Abnormal measurement data. Please do not wear shoes and socks for measurement, and remove metal accessories and other metal jewel from your body.
- Long continuous use triggers overheating protection. Please wait for the product to cool down before reusing it
- To reset to factory settings, please press and hold the left and right selection keys for about simultaneously **A3** (about 10 seconds) to reset until the “HELLO” is appearing.

VII. STORING

After using the scale, put it into the horizontal position (not vertically), out of reach of children and incapacitated people.

VIII. MAINTENANCE


Do not use coarse and aggressive cleaning agents (e.g. **sharp objects, scrapers, diluting agents or other solvents**). Clean the board and the handheld sensor by wiping with a soft wet cloth. Make sure that water does not get inside the appliance.

More extensive maintenance or maintenance requiring intervention in the inner parts of the appliance must be carried out by professional service!

Manufacturer contact: info@eta.cz

Failure to comply with the instructions of the manufacturer will lead to expiration of the right to repair!

The information regarding warranty and product repairs you can find on the website www.eta.cz/servis-eu.

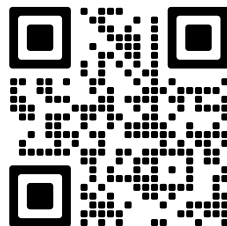IX. ENVIRONMENTAL PROTECTION

If the dimensions allow, there are labels of materials used for production of packing, components and accessories as well as their recycling on all parts. The symbols specified on the product or in the accompanying documentation mean that the used electrical or electronic products must not be disposed of together with municipal waste. For proper disposal of the product, hand it over at special collection places where it will be taken over free of charge. Suitable disposal of the product can help to maintain valuable natural resources and to prevent possible negative impacts on the environment and human health, which could be the possible consequences of improper waste disposal. Fines can be imposed for improper disposal of this type of waste in agreement with the national regulations. This appliance is equipped with Li-ion batteries with a long service life. For environmental protection it is necessary to remove the battery from the appliance after its service life and the battery and the appliance must be disposed of safely using special collection networks. Ask for more details at the local authority or at the local waste collection site.

Removing the battery

Remove the battery from the appliance only when completely discharged. Dismantle the scale. Gradually disconnect and remove the batterie.

X. TECHNICAL DATA

Weighing range max.:	180 kg / 396 lb / 28 st
Discrimination:	<100 kg 50 g / >100 kg 100 g
Net weight of device:	2.7 kg
Size of the product (mm):	375 x 352 x 35
Protection class of the appliance:	III.
Frequency range:	2402 - 2480 MHz
Output power (max):	3.4 dBm
Bluetooth version:	5.0
Compatible operating systems:	Android 4.3; iOS 9.0 (and higher)
Maximum communication reach:	10 m (open space without any obstacles and interferences).
WIFI network -	2,4 GHz WPA1 and WPA2 security (recommended WPA2)

NOTICE

- The device cannot be operated on IEEE 802.1X-secured Wifi networks (usually corporate Wifi networks).

Hereby, ETA a.s. declares that the radio equipment type ETA8779 is in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the following internet address: <http://eta.cz/declaration-of-conformity>

The manufacturer reserves the right to change the technical specifications for the respective model. The manufacturer reserves the right to change/update the ETA Body Balance application without prior notice. Due to the character of such changes/updates, some information or procedures in this manual might become obsolete or might be missing.

The capacity of the product battery decreases over time. Decreasing of battery capacity occurs during normal use of the appliance and is considered usual wear and tear of the appliance. Decreasing battery capacity is not a defect of the appliance.

WARNINGS AND SYMBOLS USED ON THE APPLIANCE, PACKAGING OR IN THE INSTRUCTIONS MANUAL:

HOUSEHOLD USE ONLY. DO NOT IMMERSE IN WATER OR OTHER LIQUIDS.

TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, ARRANGES OR PLAYPENS. THIS BAG IS NOT A TOY.

WARNING

It is not allowed to dispose it together with communal waste

Please read the instructions manual

Removable power supply unit

Do not step onto the scale with wet feet or on wet surface of the scale, there is a risk of slipping.

AUTO STEP ON – just step on the scale and stand still, the scale will automatically turn ON and start taking measurements.

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© DATE 05/06/2025
Service: Omega electric GmbH, Servicezentrale, Industriering 2, D- 04626 Schmölln.